




Pliers ■

State of art technology

**High quality pliers without plating,
perfect performance and quality guarantee.**


VST Orthodontic Pliers

VST ORTHODONTIC PLIERS

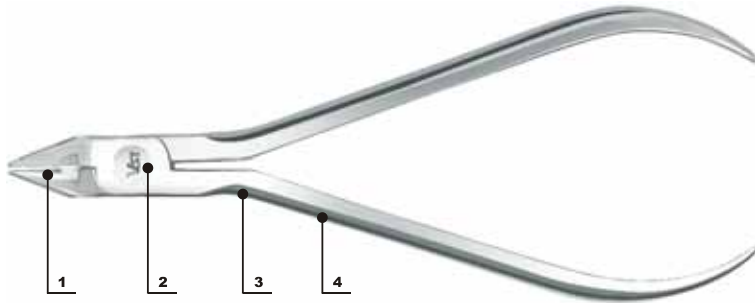

Pliers Features:

1. The new VST series has solved the problem of tip's color changing and tip's broken by importing international advanced technology.
2. Specially designed zero clearance hinge makes handles more tightly connected, and would not be loosened during operation.
3. Designed in accordance with ergonomics and rounded edges, makes dentists and patients more safety and comfortable.
4. Fine imported medical stainless steels, the pliers has been carefully grounded and polished, perfect in workmanship, excellent corrosion resistant and heat-resistant.
5. Made by CNC production lines with exquisite fixtures and moulds, ensuring precision and high quality.

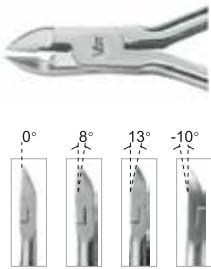

Q106-01	0°	Standard	1pc/pk
Q106-01L	0°	Long handle	1pc/pk
Q106-01M	0°	Mini	1pc/pk
Q106-02	8°	Standard	1pc/pk
Q106-02L	8°	Long handle	1pc/pk
Q106-03	13°	Standard	1pc/pk
Q106-03L	13°	Long handle	1pc/pk
Q106-04L	-10°	Long handle	1pc/pk

Light wire cutter

Hardness : HRC58-63,
Cutting edge made by new Technology,
ligature wires or elastic rings' diameter
up to 0.30mm (0.021")
(can not be used for cutting arch wires)



Q106-05	Standard	1pc/pk
Q106-05L	Long handle	1pc/pk
Q106-05M	Mini	1pc/pk
Q106-05LM	Long handle Mini	1pc/pk

Distal end cutter

Hardness : HRC 58-63,
Cutting edge made by new technology.
Used to cut the long distal end of arch
wires, and safety holding the cut part.
Arch wires up to 0.021" X0.025"
(0.53mmX0.64mm)



Q106-06	Standard	1pc/pk
Q106-06M	Mini	1pc/pk
Q106-06L	Long handle	1pc/pk
Q106-06LM	Long handle Mini	1pc/pk

**Distal end cutter
- with flush end**

Hardness : HRC 58-63,
Cutting edge made by new technology.
Used to cut the distal end of arch wires
and safety holding the cut part, the distal
end of arch wire will be flush with tube
slot' s distal end. Arch wires up to
0.021" X0.025" (0.53mmX0.64mm)

VST Orthodontic Pliers

Light wire pliers

Hardness of tip: HRC 55-60,
Cutting edge made by new Technology,
Used to make different fine loops for arch wires,
and to make spring loops. Wires' diameter up
to 0.030" (0.76mm)

Q106-07	Standard	1pc/pk
Q106-07L	Long handle	1pc/pk



Light wire pliers with cutter

Hardness of tip: HRC55-60,
Cutting edge made by new technology,
Used to make different fine loops for
arch wires, and to make spring loops.
Tip with cutting edge for wires cutting.

Q106-08	Standard	1pc/pk
Q106-08L	Long handle	1pc/pk



Heavy cutter

Hardness of tip: HRC 58-63,
Cutting edge made by new technology,
Used to cut Australia wires, NITI wires,
stainless steel wires, Wires' diameter up
to 0.030" (0.76mm)

Q106-09	Standard	1pc/pk
Q106-09L	Long handle	1pc/pk



NITI cinch back pliers

Hardness of tip: HRC55-60,
Cutting edge made by new technology.
Used to cinch back distal end of arch
wires to avoid any injury of soft tissue.
NITI wires' diameter up to 0.021" (0.53mm)

Q106-10	Dual direction	Standard	1pc/pk
Q106-10L	Dual direction	Long handle	1pc/pk
Q106-11L	One direction left	Long handle	1pc/pk
Q106-12L	One direction right	Long handle	1pc/pk



Left



Right



V stop pliers

To make V stop along arch
wires to prevent any movement.

Q107-01	1pc/pk
---------	--------



Pliers ■

VST Orthodontic Pliers



Q107-02	Standard	1pc/pk
Q302-01	Replacement tips	10pcs/pk

Cap removing pliers

To remove the caps of convertible tubes.



Q107-03	1pc/pk
---------	--------

Multi-function pliers

Hardness: HRC 50-55
Used to bend, flatten and cut wires.



Q105-01	1pc/pk
---------	--------

Torque pliers


Hardness of tip: HRC55-60 ,
Cutting edge made by new technology.
Used in pair to make torques on rectangle wires, thickness of tips: 0.050" (1.27mm)



Q105-02	1pc/pk
---------	--------

Tweed loop pliers

Hardness of tip: HRC55-60 ,
Cutting edge made by new technology.
To make small loops, distal stops and Ω loops on arch wires

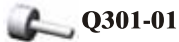

Q105-16	Posterior/short tip	1pc/pk
Q105-04	Anterior	1pc/pk
Q301-01	Replacement tips	10pcs/pk

Band remover

Hardness of tip: HRC50-55 ,
Used to remove bands and other bondable attachment.



Q105-05	1pc/pk
---------	--------

Bracket removing pliers

Hardness of tip: HRC50-55 ,
Used to remove anterior brackets.

VST Orthodontic Pliers

Posterior bracket removing pliers

Hardness of tip: HRC50-55 ,
Used to remove posterior brackets.

Q105-06	Standard	1pc/pk
Q105-06L	Long handle	1pc/pk



Weingart pliers

Hardness of tip: HRC50-55 ,
Used to cinch back distal end of arch wires
to avoid any injury of soft tissue.

Q105-07	Standard	1pc/pk
Q105-07L	Long handle	1pc/pk



Crimpable hook pliers

Hardness of tip: HRC50-55 ,
To fix the crimpable hooks along arch wires.

Q105-08	Standard	1pc/pk
Q105-08L	Long handle	1pc/pk



Loop bending pliers

Hardness of tip: HRC50-55 ,
To make different small loops for
multi-loop edgewise arch wires,
diameter of loops: 2.00, 3.20, 4.20 mm.

Q105-09	1pc/pk
---------	--------



How' s pliers

Hardness of tip: HRC50-55 ,
To hold arch wires for guiding
it into buccal tubes, to hold ligature
wires, or to form anterior band.

Q105-10	1pc/pk
---------	--------



Kim pliers

Hardness of tip: HRC50-55 ,
To bend multi-loop edgewise
wires, with cutter on tips.

Q105-11	1pc/pk
---------	--------



Pliers ■

VST Orthodontic Pliers



Q105-12	1pc/pk
----------------	--------

Nance loop bending pliers

Hardness of tip: HRC50-55 ,
To make vertical loops of different height along arch wires, 3.0, 4.0, 5.0, 6.0mm
Wires' diameter up to : 0.022" (0.56mm)



Q105-13	Standard	1pc/pk
Q105-13L	Long handle	1pc/pk
Q301-02	Replacement tips	10pcs/pk
Q302-02	Replacement blade	10pcs/pk

Bondage removing pliers

Hardness of tip: HRC50-55 ,
Used to remove bondage on the crown without any injury on enamel.



Q301-02 Q302-02



Q105-14	1pc/pk
----------------	--------

Three jaw pliers

Hardness of tip: HRC50-55 ,
To make V loops and bend arch wires.
Wires' diameter up to 0.030" (0.76mm)



Q105-15	1pc/pk
----------------	--------

Omega loop bending pliers

Hardness of tip: HRC50-55 ,
To make fine Ω loops. Three different loops :
diameter 0.063" , 0.075" , 0.10"
(diameter 1.6mm, 2.0mm, 2.6mm)



Q103-01	1pc/pk
----------------	--------

Elastic plier

Hardness of tip: HRC50-55 ,
For expanding elastic separating rings before separating molar gap.

VST Orthodontic Pliers

Packaging box

Including guarantee sheet
 Designed to pack pliers to avoid any damage during shipment and storage.



Pliers stands

Clear, modified design
 Stainless steel bar, high temperature sterilization not permitted. plastic lateral plates.

Q201-03		1set/pk
----------------	--	---------



Fully stainless steel stands
 Newly designed, high temperature sterilized.

Q201-04	Standard	1set/pk
Q201-05	Large	1set/pk

